

## REF 76191212

spitzenlos

4-Draht-Körbchen, gerade

tipless


4-wire-basket, straight

### Der »STONIZER®« advanced

Mit kreisrunder Korbform und neuer Schlaufentechnik.

### The »STONIZER®« advanced

With circular basket design and new loop technique.

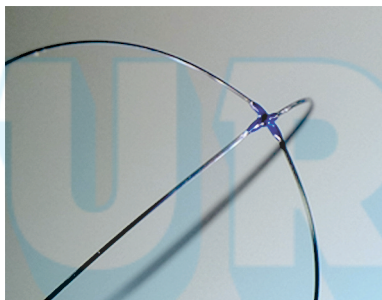

### Fangkörbchen zur ureterorenoskopischen Steinentfernung

- zur Anwendung in der Niere geeignet
- knickstabil
- bequeme Einhand-Bedienung
- kreisrunde Korbform erleichtert Einfangen des Steines
- neuartige Konstruktion des Korbes und modifizierte Materialien sichern höchste Zuverlässigkeit – auch beim Bergen von Steinstraßen
- enorme Zug- und Druckstabilität des Körbchens durch beständiges Material
- Mechanismus zur Begrenzung der Zugkraft gewährleistet sichere Steinentfernung
- verbesserte Sicht durch reflektieren des Korbende
- atraumatisches distales Korbende
- optimale Schließfunktion
- Steuerung des Körbchens mit einer hohen Drehstabilität
- hohe Funktionalität auch bei abgewinkeltem Schaft
- flexibel durch Nitinol-Material
- erhöhter Spülfluss im Ureterorenoskop durch geringen Durchmesser
- abnehmbarer und wieder montierbarer Handgriff durch Rändelschraube
- zur einmaligen Verwendung
- Einführhilfe

### Retrieval basket for ureterorenoscopic stone extraction

- for use in kidneys
- kink resistant
- comfortable one-handed-use
- circular basket design facilitates capture of stones
- innovative construction of basket and modified materials ensure the highest possible reliability – also in case of stone streets
- enormous pull and pressure stability of the basket because of its durable material
- mechanism to limit the tensile strength provides secure stone retrieval
- improved visibility by reflecting end of basket
- atraumatic distal end of basket
- optimal closing function
- guidance of the basket with high torque stability
- high functionality even with bended shaft
- flexible Nitinol material
- increased irrigation in the ureterorenoscope due to small diameter
- removable and re-mountable handle by knurled screw
- single use
- insertion guide

REF REF	Größe (FR/CH) Size (FR/CH)	Länge Length	Ø Korb Ø Basket
76191212	Ch. 1,9	120 cm	12 mm

Verpackung: steril  
Original-Karton: 1 Stück

Packing: sterile  
Original carton: 1 piece