

# Bendable


## Ureteral Access Suction Sheath


### Flexible Bendable Tip

Easier Search And Closer to Renal Calyx

### Controllable Suction

Safety Intraluminal Pressure With Higher Stone Free Rate

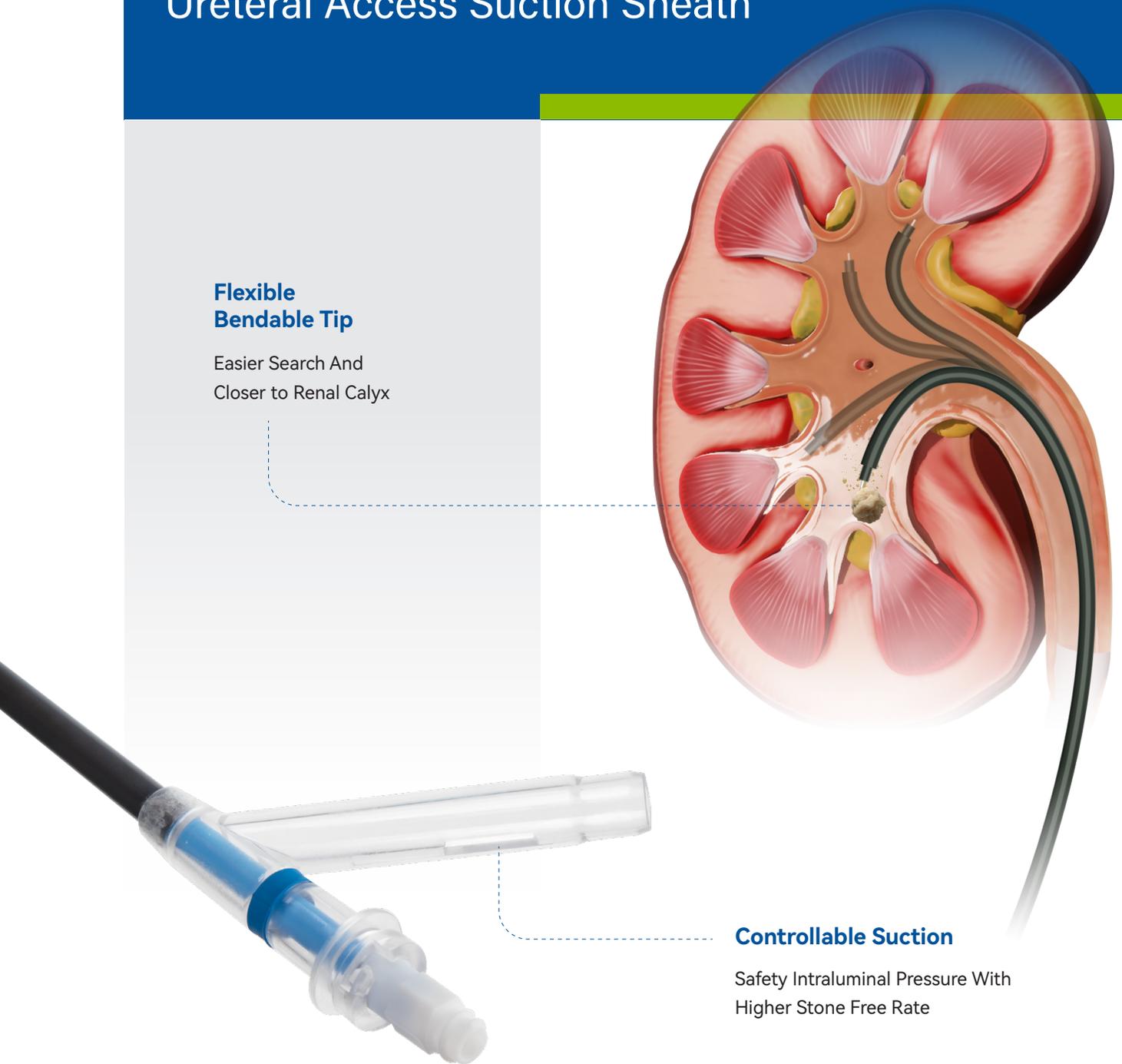

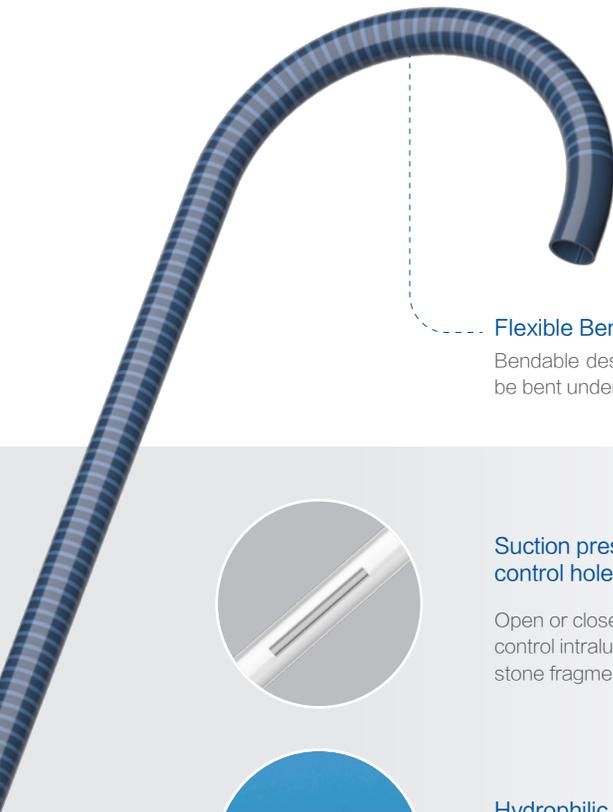

### Flexible Bendable Tip

Bendable design distal end of sheath can be bent under the control of FURS.



### Suction pressure control hole

Open or close the lateral suction hole to control intraluminal pressure and suck out stone fragments



### Hydrophilic coating

Allows for ease of insertion. Improved coating is designed for durability in bilateral cases.




### Stone filtration And Collection

A filter is designed to collect fragments and prevent suction tube from being blocked



### Core Reinforced

The core consists of a patented reinforced coil construction to provide optimal flexibility and maximum resistance to kinking and compression.

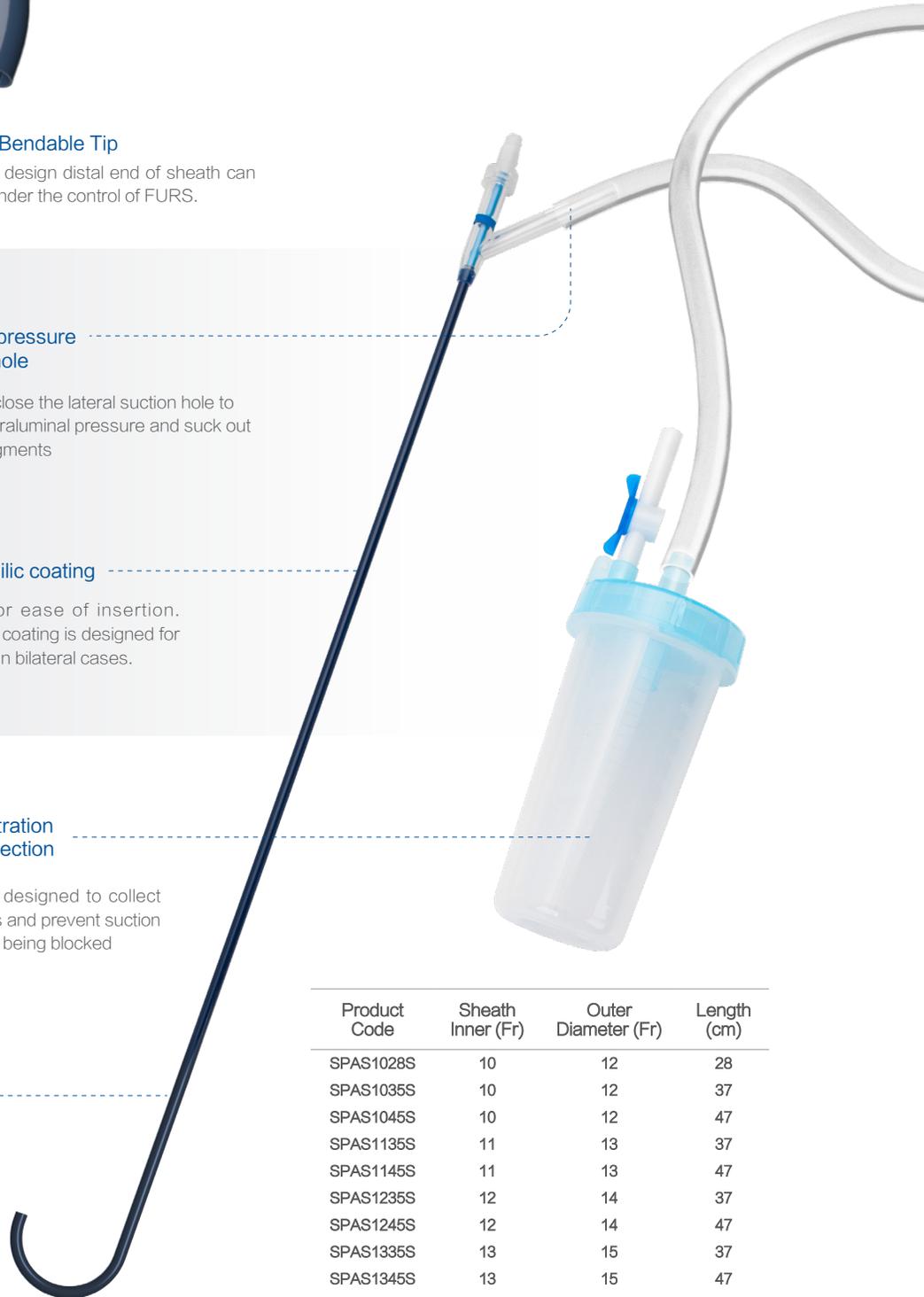

Product Code	Sheath Inner (Fr)	Outer Diameter (Fr)	Length (cm)
SPAS1028S	10	12	28
SPAS1035S	10	12	37
SPAS1045S	10	12	47
SPAS1135S	11	13	37
SPAS1145S	11	13	47
SPAS1235S	12	14	37
SPAS1245S	12	14	47
SPAS1335S	13	15	37
SPAS1345S	13	15	47
SPAS1435S	14	16	37
SPAS1445S	14	16	47