


Guidewire

Zebra Nitinol Core Guidewire

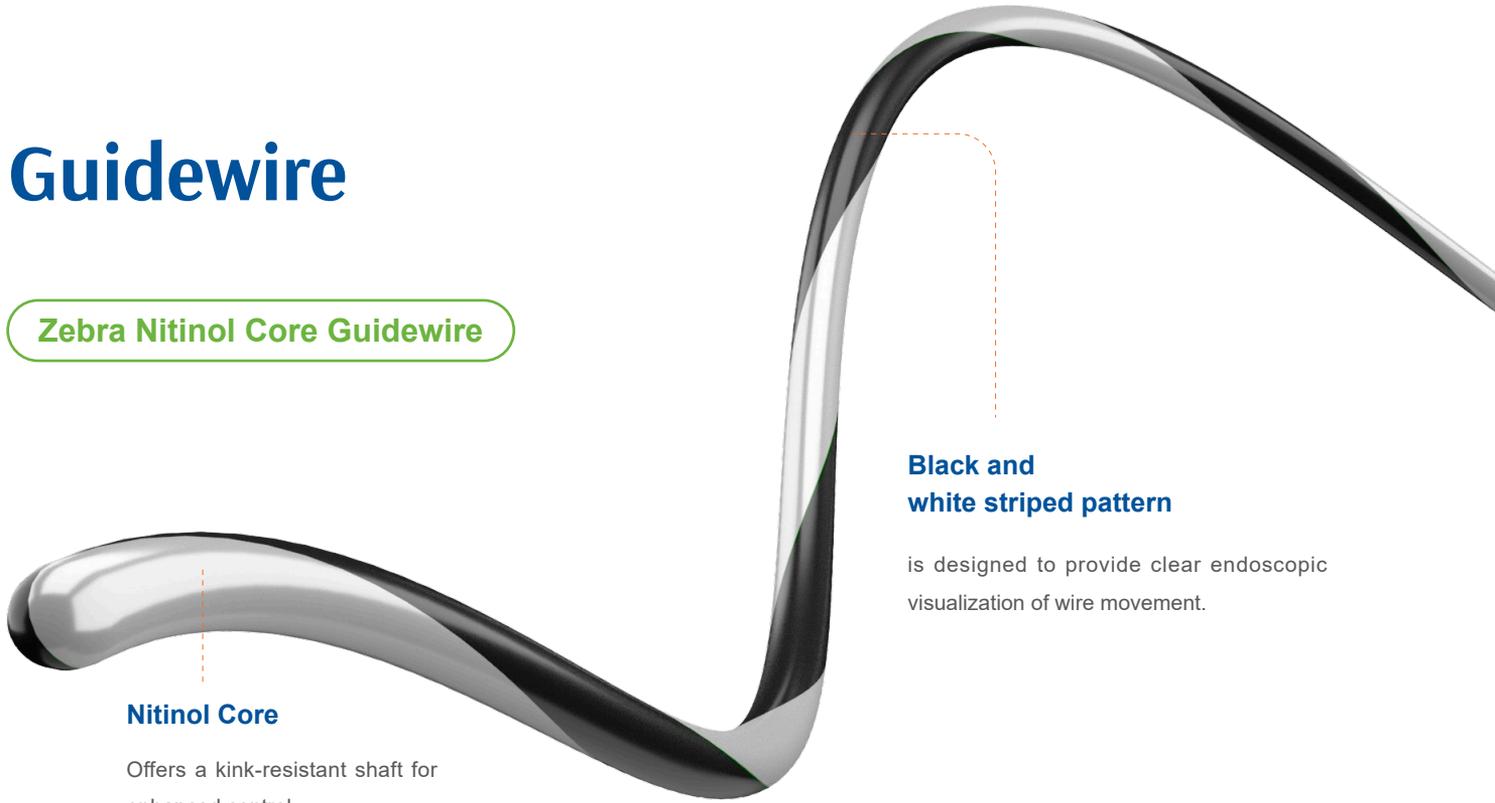

Nitinol Core

Offers a kink-resistant shaft for enhanced control.

Black and white striped pattern

is designed to provide clear endoscopic visualization of wire movement.

Multiple Tip Options

Designed to provide the additional tip flexibility for negotiation of narrowed ureters and access beyond obstructions.

Product Code	Length (cm)	Size	Tip
SPGW3232	150	ø0.032"	Angled
SPGW3235	150	ø0.035"	Angled
SPGW3132	150	ø0.032"	Straight
SPGW3135	150	ø0.035"	Straight

Insertor

- 01 SCOOP Proximal end has a large open trough to allow easy insertion of wire guides and other instrumentation.
- 02 GRIP Exterior grooves provide a convenient grip for additional control and precise insertion.
- 03 TIP Distal tip has a slender taper to ensure a snug fit. Can also be used to puncture the seal on an endoscopic valve. Insertor is compatible with various adapters.

Product Code	Accepts Instrumentation (Fr)	Length(cm)
SPIS1100	3.1	5

